
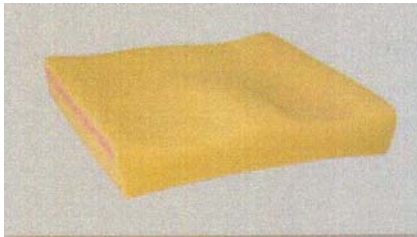


# PRISM Seating & Positioning


## Ideal Cushion

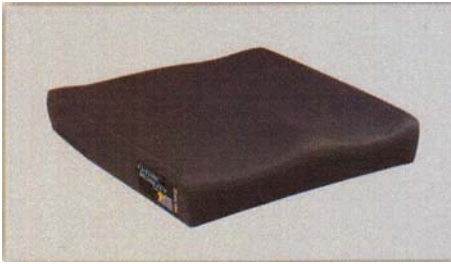

The Prism Ideal Cushion was designed for users that require mild to moderate postural support and offers gentle contours helping to promote correct seating posture.

A seat cushion with 3 distinct layers: a top layer of soft, comfortable poly foam; a central layer of pressure relieving visco-elastic foam; and a base layer of high density polyurethane support foam contour cut to provide positioning.

We have created an economical and lightweight cushion with comfort, positioning and pressure management. These qualities make the Ideal Cushion just that.

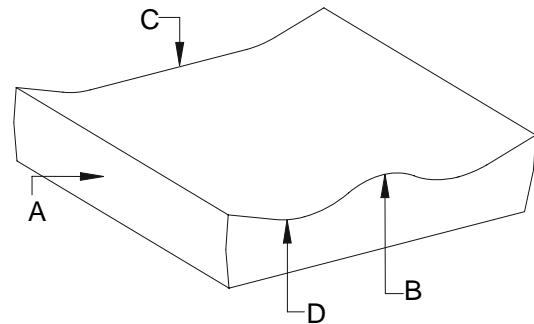

### Prism™ Ideal Cushion

#### Stock Sizes—Next Day Shipments

16" x 16"	IC1616	16" x 18"	IC1618
18" x 16"	IC1816	18" x 18"	IC1818
20" x 16"	IC2016	20" x 18"	IC2018

#### Non Stock Sizes—Shipments within 3 business days

16" x 20"	IC1620	17" x 16"	IC1716
17" x 18"	IC1718	17" x 20"	IC1720
18" x 20"	IC1820	20" x 20"	IC2020

### CUSHION THICKNESS

- A—Height at adductor = 3"
- B—Height at abductor = 3"
- C—Height at rear IT = 2 1/4"
- D—Leg trough = 2 1/4"

### SPOTLIGHT ON

- Waterfall front edge for foot propelling comfort
- Low profile high-density foam base enables proper positioning
- Core layer of visco-elastic foam to provide comfort and pressure reduction
- Light weight cushion

### WARRANTY

Two (2) years limited on Cushion and 90 days on cover.

### PRODUCT WEIGHT

2 lbs

### USER WEIGHT

250 lbs.

