
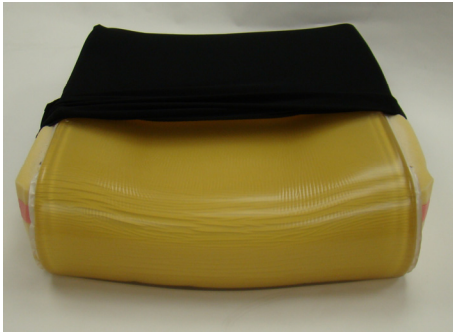


New

PRISM Seating & Positioning

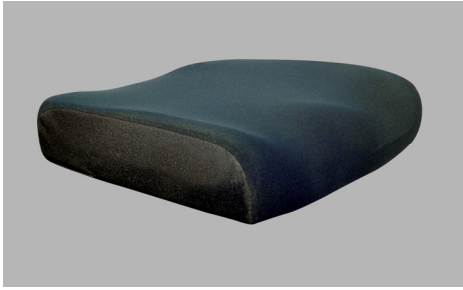


Supreme Gel Cushion

The Prism Supreme Gel Cushion was designed for users requiring mild to moderate positioning support and excellent pressure distribution relief.

The gel cube overlay reduces shear and friction that contribute to skin breakdown, while the anatomically designed firm polyurethane base provides superior pelvic stability and leg positioning.

A central core of pressure relieving visco-elastic insert under the IT's allows for excellent immersion for bony prominences.

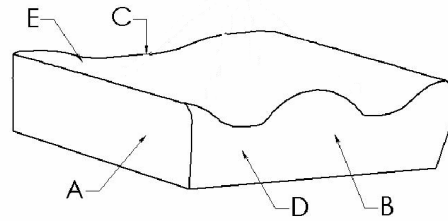

SPOTLIGHT ON

- Gel overlay provides even distribution of weight across the cushion
- Waterfall front edge for foot propelling comfort
- Cane cut in rear to allow for correct placement in wheelchairs and proper seat depth
- Radius cut to match shape of back upholstery
- 3" of visco foam under IT's for superior immersion
- Supreme pressure relieving
- Standard water resistant inner cover
- Choose from our standard high 4-way stretch outer cover

CUSHION THICKNESS

- A—Height at adductor = 5"
- B—Height at abductor = 5"
- C—Height at rear IT = 3 1/4"
- D—Leg trough = 3 1/4"
- E— Height at rear trochanter = 4 1/4"

WARRANTY

Two (2) years limited on Cushion and 90 days on cover.

PRODUCT WEIGHT

10 lbs

USER WEIGHT

250lbs

Prism Supreme Cushion

Stock Sizes—same day shipments

14" x 14"	SGC1414	14" x 16"	SGC1416
16" x 16"	SGC1616	16" x 18"	SGC1618
18" x 16"	SGC1816	18" x 18"	SGC1818
18" x 20"	SGC1820		

Non-stocked Sizes—shipments within 3 business days

15" x 16"	SGC1516	15" x 18"	SGC1518
16" x 20"	SGC1620	17" x 16"	SGC1716
17" x 18"	SGC1718	17" x 20"	SGC1720
19" x 16"	SGC1916	19" x 18"	SGC1918
19" x 20"	SGC1920	20" x 16"	SGC2016
20" x 18"	SGC2018	20" x 20"	SGC2020



3223 Orlando Drive, Mississauga, ON, L4V 1C5, 905-671-1661, 888-737-4011, fax 866-670-3377 www.future-mobility.com
Future Mobility Healthcare Inc. Orion and Prism are trademarks of Future Mobility Healthcare Inc.