

New



Future Mobility Healthcare Inc., manufactures and markets products that increase the freedom and safety of mobility challenged individuals. Reliable for quality and service, Future Mobility Healthcare Inc. provides high quality and innovative mobility products to the home healthcare market.



SPOTLIGHT ON:

- Easy to use 2 button hand control with flexible cable
- Alarm sounds when batteries are low and require recharging
- Supports up to 250 lbs occupant


FEATURES:

- 3 to 45 degrees of infinite tilt
- Available on Orion II chairs with seat width of 15" to 20"
- Seat-to-floor heights: 13 1/4" - 19"
- 27V Plug in Battery Charger (120 V AC)
- 24V High Capacity Battery/Control Box on Board
- Compact Actuator with 135mm (5.31") Stroke Length



**24V High Capacity
Battery/Control Box**

Compact Actuator

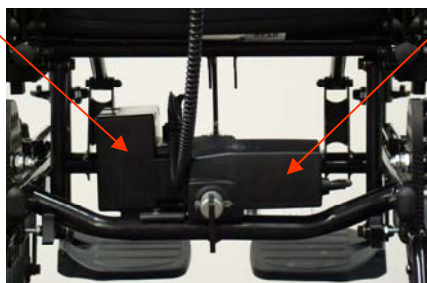

Hand Controller - HB40

- High strength plastic housing
- Ergonomic design & easy touch buttons
- Hook for easy placement

WARRANTY

Power Tilt batteries are warranted for six(6) months while the control box, battery charger, actuator and hand control are warranted for one (1) year.