

Supreme Cushion



The Prism Supreme Cushion was designed for users requiring moderate to high positioning support and provides excellent pressure distribution relief.

The anatomically designed firm polyurethane base with its aggressive leg troughs and abductor contouring provides superior pelvic stability and leg positioning. A central core of pressure relieving visco-elastic insert under the IT's allows for excellent immersion for bony prominences.

Supreme comfort, supreme positioning and supreme pressure relief—that's the Prism Supreme!

SPOTLIGHT ON

- High density, contour cut layered foam for a consistent, long lasting product
- Waterfall front edge for foot propelling comfort
- Cane cut in rear to allow for correct placement in wheelchairs and proper seat depth
- Radius cut to match shape of back upholstery
- 3" of visco foam under IT's for superior immersion
- Pressure relieving
- Standard water resistant inner cover
- Choose from our standard high 4-way stretch outer cover

WARRANTY


Two (2) years limited on Cushion and 90 days on cover.

PRODUCT WEIGHT

3.45lbs

USER WEIGHT

250lbs

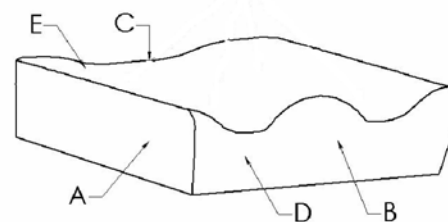

CUSHION THICKNESS

- A—Height at adductor = 5"
 - B—Height at abductor = 5"
 - C—Height at rear IT = 3 1/4"
 - D—Leg trough = 3 1/4"
 - E— Height at rear trochanter = 4 1/4"
- Air inserts = adjustable areas

Prism Supreme Cushion

Stock Sizes—same day shipments

14" x 14"	SC1414	14" x 16"	SC1416
16" x 16"	SC1616	16" x 18"	SC1618
18" x 16"	SC1816	18" x 18"	SC1818
18" x 20"	SC1820		

Non-stocked Sizes—shipments within 3 business days

15" x 16"	SC1516	15" x 18"	SC1518
16" x 20"	SC1620	17" x 16"	SC1716
17" x 18"	SC1718	17" x 20"	SC1720
19" x 16"	SC1916	19" x 18"	SC1918
19" x 20"	SC1920	20" x 16"	SC2016
20" x 18"	SC2018	20" x 20"	SC2020

